
- en **Instructions for use/Technical description**
Laryngeal mirror and laryngeal mirror handle
- USA Note for U.S. users
This Instructions for Use is NOT intended for United States users. Please discard. The Instructions for Use for United States users can be obtained by visiting our website at www.aesculapusa.com. If you wish to obtain a paper copy of the Instructions for Use, you may request one by contacting your local Aesculap representative or Aesculap's customer service at 1-800-282-9000. A paper copy will be provided to you upon request at no additional cost.
- de **Gebrauchsanweisung/Technische Beschreibung**
Kehlkopfspiegel und Kehlkopfspiegelgriffe
- fr **Mode d'emploi/Description technique**
Laryngoscopes et manches de laryngoscope
- es **Instrucciones de manejo/Descripción técnica**
Espejo laringeo y mango del espejo laringeo
- it **Istruzioni per l'uso/Descrizione tecnica**
Specchietto laringeo e impugnatura per specchietto laringeo
- pt **Instruções de utilização/Descrição técnica**
Espelho laringeo e punho para espelho laringeo
- nl **Gebruiksaanwijzing/Technische beschrijving**
Larynxspiegel en larynxspiegelhandgreep
- da **Brugsanvisning/Teknisk beskrivelse**
Larynxspejl og håndtag til larynxspejl
- sv **Bruksanvisning/Teknisk beskrivning**
Larynxspegel och handtag till larynxspegel
- fi **Käyttöohje/Tekninen kuvaus**
Laryngeaalinen peili ja laryngeaalisen peilin kahva
- lv **Lietošanas instrukcijas/tehniskais apraksts**
Balseņu spogulis un balseņu spoguļa rokturis
- lt **Naudojimo instrukcija/techninis aprašas**
Gerklų veidrodėlis ir gerklų veidrodėlio rankena
- ru **Инструкция по применению/Техническое описание**
Гортанное зеркало и рукоятка для гортанного зеркала
- cs **Návod k použití/Technický popis**
Laryngoskopické zrcátko a rukojeť laryngoskopického zrcátka
- pl **Instrukcja użytkowania/Opis techniczny**
Lusterko krtańowe i uchwyty do lusterka krtańowego
- sk **Návod na použitie/Technický opis**
Hrtanové zrkadlo a rúčka hrtanového zrkadla
- hu **Használati útmutató/Műszaki leírás**
Gége-tükör és gégetükörnyél
- sl **Navodila za uporabo/Tehnični opis**
Laringealno ogledalce in ročaj laringealnega ogledalca
- hr **Upute za uporabu/Tehnički opis**
Laringealno ogledalo i ručka laringealnog ogledala
- ro **Manual de utilizare/Descriere tehnică**
Oglindă laringiană și mâner pentru oglindă laringiană
- bg **Упътване за употреба/Техническо описание**
Ларингеално огледало и ръкохватка за ларингеално огледало
- tr **Kullanım Kılavuzu/Teknik açıklama**
Larenks aynası ve larenks aynası sapı
- el **Οδηγίες χρήσης/Τεχνική περιγραφή**
Λαρυγγικό κάτοπτρο και λαβή λαρυγγικού κατόπτρου

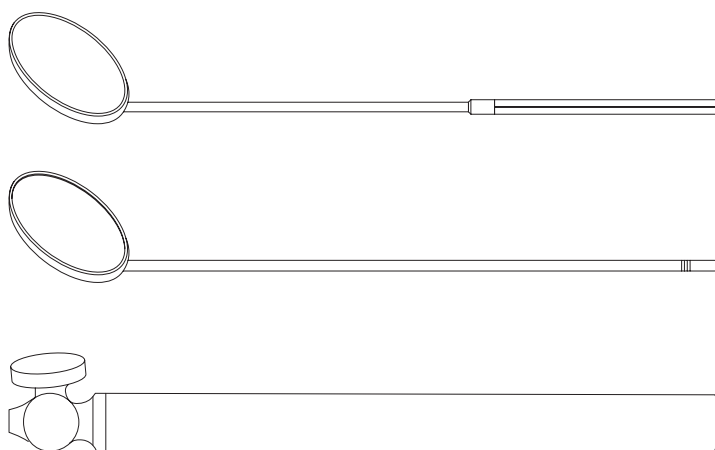1. About this document

Note

General risk factors associated with surgical procedures are not described in these instructions for use.

1.1 Scope

These instructions for use apply for the following products:

1.1.1 Laryngeal mirror with handle

Art. no.	Designation
OM078R	Hexagonal laryngeal mirror with handle Ø 6 mm
OM079R	Hexagonal laryngeal mirror with handle Ø 8 mm
OM080R	Hexagonal laryngeal mirror with handle Ø 10 mm
OM081R	Hexagonal laryngeal mirror with handle Ø 12 mm
OM082R	Hexagonal laryngeal mirror with handle Ø 14 mm
OM083R	Hexagonal laryngeal mirror with handle Ø 16 mm
OM084R	Hexagonal laryngeal mirror with handle Ø 18 mm
OM085R	Hexagonal laryngeal mirror with handle Ø 20 mm
OM086R	Hexagonal laryngeal mirror with handle Ø 22 mm
OM087R	Hexagonal laryngeal mirror with handle Ø 24 mm
OM088R	Hexagonal laryngeal mirror with handle Ø 26 mm
OM089R	Hexagonal laryngeal mirror with handle Ø 28 mm
OM090R	Hexagonal laryngeal mirror with handle Ø 30 mm

1.1.2 Laryngeal mirror without handle

Art. no.	Designation
AC468R	Laryngeal mirror Ø 6 mm
AC469R	Laryngeal mirror Ø 8 mm
AC470R	Laryngeal mirror Ø 10 mm
AC471R	Laryngeal mirror Ø 12 mm
AC472R	Laryngeal mirror Ø 14 mm
AC473R	Laryngeal mirror Ø 16 mm
AC474R	Laryngeal mirror Ø 18 mm
AC475R	Laryngeal mirror Ø 20 mm
AC476R	Laryngeal mirror Ø 22 mm
AC477R	Laryngeal mirror Ø 24 mm
AC478R	Laryngeal mirror Ø 26 mm
AC479R	Laryngeal mirror Ø 28 mm
AC480R	Laryngeal mirror Ø 30 mm

1.1.3 Handle for laryngeal mirror

Art. no.	Designation
AC481R	Handle for laryngeal mirror 110 mm

► For article-specific instructions for use as well as information on material compatibility and lifetime see B. Braun eIFU at eifu.bbraun.com

2. Safety messages

Safety messages make clear the dangers to patient, user and/or product that could arise during the use of the product. Safety messages are labeled as follows:

⚠ WARNING

Indicates a possible threat of danger. If not avoided, minor or moderate injury may result.

⚠ CAUTION

Indicates a possible threat of material damage. If not avoided, the product may be damaged.

2. Clinical use

2.1 Areas of use and limitations of use

2.1.1 Intended use

Laryngeal mirror grips serve to receive interchangeable laryngeal mirrors. This mirror is used for indirect viewing of the larynx.

2.1.2 Indications

Note

The manufacturer is not responsible for any use of the product against the specified indications and/or the described applications.

For indications, see Intended use.

2.1.3 Contraindications

No known contraindications.

2.2 Safety information

2.2.1 Clinical user

General safety information

To prevent damage caused by improper setup or operation, and to not compromise the manufacturer warranty and liability:

- Use the product only according to these instructions for use.
- Follow the safety and maintenance instructions.
- Ensure that the product and its accessories are operated and used only by persons with the requisite training, knowledge and experience.
- Store any new or unused products in a dry, clean, and safe place.
- Prior to use, check that the product is in good working order.
- Keep the instructions for use accessible for the user.

Note

The user is obligated to report all severe events in connection with the product to the manufacturer and the responsible authorities of the state in which the user is located.

Notes on surgical procedures

It is the user's responsibility to ensure that the surgical procedure is performed correctly.

Appropriate clinical training as well as a theoretical and practical proficiency of all the required operating techniques, including the use of this product, are prerequisites for the successful use of this product.

The user is required to obtain information from the manufacturer if there is an unclear preoperative situation regarding the use of the product.

2.2.2 Sterility

The product is delivered in an unsterile condition.

- Clean the new product after removing its transport packaging and prior to its initial sterilization.

2.3 Application

⚠ WARNING

Risk of injury and/or malfunction!

- Prior to each use, inspect the product for loose, bent, broken, cracked, worn, or fractured components.
- Always carry out a function test prior to each use of the product.

3. Validated reprocessing procedure

3.1 General safety instructions

Note

Adhere to national statutory regulations, national and international standards and directives, and local, clinical hygiene instructions for sterile processing.

Note

For patients with Creutzfeldt-Jakob disease (CJD), suspected CJD or possible variants of CJD, observe the relevant national regulations concerning the reprocessing of products.

Note

Mechanical reprocessing should be favored over manual cleaning as it gives better and more reliable results.

Note

Successful processing of this medical device can only be ensured if the processing method is first validated. The operator/sterile processing technician is responsible for this.

Note

If there is no final sterilization, then a virucidal disinfectant must be used.

Note

For up-to-date information about reprocessing and material compatibility, see B. Braun eIFU at eifu.bbraun.com. The validated steam sterilization procedure was carried out in the Aesculap sterile container system.

3.2 General information

Dried or affixed surgical residues can make cleaning more difficult or ineffective and lead to corrosion. Therefore the time interval between application and processing should not exceed 6 h; also, neither fixating pre-cleaning temperatures >45 °C nor fixating disinfecting agents (active ingredient: aldehydes/alcohols) should be used.

Excessive measures of neutralizing agents or basic cleaners may result in a chemical attack and/or to fading and the laser marking becoming unreadable visually or by machine for stainless steel.

Residues containing chlorine or chlorides e.g. in surgical residues, medicines, saline solutions and in the service water used for cleaning, disinfection and sterilization will cause corrosion damage (pitting, stress corrosion) and result in the destruction of stainless steel products. These must be removed by rinsing thoroughly with demineralized water and then drying.

Additional drying, if necessary.

Only process chemicals that have been tested and approved (e.g. VAH or FDA approval or CE mark) and which are compatible with the product's materials according to the chemical manufacturers' recommendations may be used for processing the product. All the chemical manufacturer's application specifications must be strictly observed. Failure to do so can result in the following problems:

- Optical changes of materials, e.g. fading or discoloration of titanium or aluminum. For aluminum, the application/process solution only needs to be of pH >8 to cause visible surface changes.
- Material damage such as corrosion, cracks, fracturing, premature aging or swelling.
- Do not use metal cleaning brushes or other abrasives that would damage the product surfaces and could cause corrosion.
- Further detailed advice on hygienically safe and material-/value-preserving reprocessing can be found at www.ak-i.org, link to "AKI-Brochures", "Red brochure".

3.3 Reusable products

Influences of the reprocessing which lead to damage to the product are not known.

A careful visual and functional inspection before the next use is the best opportunity to recognize a product that is no longer functional, see Inspection.

3.4 Preparations at the place of use

- If applicable, rinse non-visible surfaces preferably with deionized water, with a disposable syringe for example.
- Remove any visible surgical residues to the extent possible with a damp, lint-free cloth.
- Transport the dry product in a sealed waste container for cleaning and disinfection within 6 hours.

3.5 Preparing for cleaning

- AC481R: Disassemble the handle from the mirror prior to cleaning.

3.6 Cleaning/Disinfection

3.6.1 Product-specific safety information on the reprocessing method

Damage to or destruction of the product due to inappropriate cleaning/disinfecting agents and/or excessive temperatures!

- Following the manufacturer's instructions, use cleaning and disinfecting agents
 - that are approved for plastic material and high-grade steel,
 - that do not attack softeners (e.g., in silicone).
- Observe specifications regarding concentration, temperature and exposure time.
- Do not exceed the maximum allowable disinfection temperature of 95 °C.

3.6.2 Validated cleaning and disinfection procedure

Validated procedure	Specific requirements	Reference
Manual cleaning with immersion disinfection ■ All products	<ul style="list-style-type: none"> ■ Suitable cleaning brush ■ Disposable syringe 20 ml ■ Drying phase: Use a lint-free cloth or medical compressed air 	Chapter Manual cleaning/disinfection and subsection: ■ Chapter Manual cleaning with immersion disinfection
Mechanical alkaline cleaning and thermal disinfection ■ All products except AC481R	<ul style="list-style-type: none"> ■ Place the product on a tray that is suitable for cleaning (avoid rinsing blind spots). 	Chapter Mechanical cleaning/disinfection and subsection: ■ Chapter Mechanical alkaline cleaning and thermal disinfecting
Manual pre-cleaning with brush and subsequent mechanical alkaline cleaning and thermal disinfection ■ AC481R	<ul style="list-style-type: none"> ■ Suitable cleaning brush ■ Disposable syringe 20 ml ■ Place the product on a tray that is suitable for cleaning (avoid rinsing blind spots). 	Chapter Mechanical cleaning/disinfection with manual pre-cleaning and subsection: ■ Chapter Manual pre-cleaning with a brush ■ Chapter Mechanical alkaline cleaning and thermal disinfecting

3.7 Manual cleaning/disinfection

- ▶ Prior to manual disinfecting, allow water to drip off for a sufficient length of time to prevent dilution of the disinfecting solution.
- ▶ After manual cleaning/disinfection, check visible surfaces visually for residues.
- ▶ Repeat the cleaning/disinfection process if necessary.

3.7.1 Manual cleaning with immersion disinfection

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemical
I	Disinfecting cleaning	RT (cold)	>15	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
II	Intermediate rinse	RT (cold)	1	-	D-W	-
III	Disinfection	RT (cold)	5	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
IV	Final rinse	RT (cold)	1	-	FD-W	-
V	Drying	RT	-	-	-	-

D-W: Drinking water
FD-W: Fully desalinated water (demineralized, microbiological, at least of drinking water quality)
RT: Room temperature
*Recommended: BBraun Stabimed fresh

- ▶ Note the information on appropriate cleaning brushes and disposable syringes, see Validated cleaning and disinfection procedure.

Phase I

- ▶ Fully immerse the product in the cleaning/disinfectant for at least 15 min. Ensure that all accessible surfaces are moistened.
- ▶ Clean the product with a suitable cleaning brush in the solution until all discernible residues have been removed from the surface.
- ▶ If applicable, brush through non-visible surfaces with an appropriate cleaning brush for at least 1 min.
- ▶ Mobilize non-rigid components, such as set screws, links, etc. during cleaning.
- ▶ Thoroughly rinse through these components with the cleaning disinfectant solution (at least five times), using a disposable syringe.

Phase II

- ▶ Rinse/flush the product thoroughly (all accessible surfaces) under running water.
- ▶ Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.
- ▶ Drain any remaining water fully.

Phase III

- ▶ Fully immerse the product in the disinfectant solution.
- ▶ Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.
- ▶ Rinse lumens at least 5 times at the beginning of the exposure time using an appropriate disposable syringe. Ensure that all accessible surfaces are moistened.

Phase IV

- ▶ Rinse/flush the product thoroughly (all accessible surfaces).
- ▶ Mobilize non-rigid components, such as set screws, joints, etc. during final rinse.
- ▶ Rinse lumens with an appropriate disposable syringe at least five times.
- ▶ Drain any remaining water fully.

Phase V

- ▶ Dry the product in the drying phase with suitable equipment (e.g. cloth, compressed air), see Validated cleaning and disinfection procedure.

3.8 Mechanical cleaning/disinfection

Note
The cleaning and disinfection device must be of tested and approved effectiveness (e.g. FDA approval or CE mark according to DIN EN ISO 15883).

Note
The cleaning and disinfection device used for processing must be serviced and checked at regular intervals.

3.8.1 Mechanical alkaline cleaning and thermal disinfecting

Machine type: single-chamber cleaning/disinfecting machine without ultrasound

Phase	Step	D [°C/°F]	t [min]	Water quality	Chemical/Note
I	Prerinse	<25/77	3	D-W	-
II	Cleaning	55/131	10	FD-W	<ul style="list-style-type: none"> ■ Concentrate, alkaline: <ul style="list-style-type: none"> - pH ~ 13 - <5 % anionic surfactant ■ working solution 0.5% <ul style="list-style-type: none"> - pH = 11*
III	Intermediate rinse	>10/50	1	FD-W	-
IV	Thermal disinfecting	90/194	5	FD-W	-
V	Drying	-	-	-	According to the program for cleaning and disinfection device

DW: Drinking water
FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

*Recommended: BBraun Helimatic Cleaner alkaline

- ▶ Check visible surfaces for residues after mechanical cleaning/disinfecting.

3.9 Mechanical cleaning/disinfection with manual pre-cleaning

Note

The cleaning and disinfection device must be of tested and approved effectiveness (e.g. FDA approval or CE mark according to DIN EN ISO 15883).

Note

The cleaning and disinfection device used for processing must be serviced and checked at regular intervals.

3.9.1 Manual pre-cleaning with a brush

Phase	Step	T [°C/°F]	t [min]	Conc. [%]	Water quality	Chemical
I	Disinfecting cleaning	RT (cold)	>15	2	D-W	Aldehyde-free, phenol-free, and QUAT-free concentrate, pH ~ 9*
II	Rinsing	RT (cold)	1	-	D-W	-

D-W: Drinking water
RT: Room temperature
*Recommended: BBraun Stabimed fresh

- ▶ Note the information on appropriate cleaning brushes and disposable syringes, see Validated cleaning and disinfection procedure.

Phase I

- ▶ Fully immerse the product in the cleaning/disinfectant for at least 15 min. Ensure that all accessible surfaces are moistened.
- ▶ Clean the product with a suitable cleaning brush in the solution until all discernible residues have been removed from the surface.
- ▶ If applicable, brush through non-visible surfaces with an appropriate cleaning brush for at least 1 min.
- ▶ Mobilize non-rigid components, such as set screws, links, etc. during cleaning.
- ▶ Thoroughly rinse through these components with the cleaning disinfectant solution (at least five times), using a disposable syringe.

Phase II

- ▶ Rinse/flush the product thoroughly (all accessible surfaces) under running water.
- ▶ Mobilize non-rigid components, such as set screws, joints, etc. during rinsing.

3.9.2 Mechanical alkaline cleaning and thermal disinfecting

Machine type: single-chamber cleaning/disinfection device without ultrasound

Phase	Step	T [°C/°F]	t [min]	Water quality	Chemical
I	Prerinse	<25/77	3	D-W	-
II	Cleaning	55/131	10	FD-W	<ul style="list-style-type: none"> ■ Concentrate, alkaline: <ul style="list-style-type: none"> - pH ~ 13 - <5 % anionic surfactant ■ working solution 0.5% <ul style="list-style-type: none"> - pH = 11*
III	Intermediate rinse	>10/50	1	FD-W	-
IV	Thermal disinfecting	90/194	5	FD-W	-
V	Drying	-	-	-	According to the program for cleaning and disinfection device

D-W: Drinking water
FD-W: Fully desalinated water (demineralized, low microbiological contamination: drinking water quality at least)

*Recommended: BBraun Helimatic Cleaner alkaline

- ▶ Check visible surfaces for residues after mechanical cleaning/disinfecting.

3.10 Inspection

- ▶ Allow the product to cool down to room temperature.
- ▶ Dry the product if it is wet or damp.

3.10.1 Visual inspection

- ▶ Ensure that all soiling has been removed. In particular, pay attention to mating surfaces, hinges, shafts, recessed areas, drill grooves and the sides of the teeth on rasps.
- ▶ If the product is dirty: repeat the cleaning and disinfection process.
- ▶ Check the product for damage, e.g. insulation or corroded, loose, bent, broken, cracked, worn or severely scratched and fractured components.
- ▶ Check the product for missing or faded labels.
- ▶ Check the products with long, slim shapes (in particular rotating instruments) for deformities.
- ▶ Check the product for damage to the spiral element.
- ▶ Check the surfaces for rough spots.
- ▶ Check the product for burrs that could damage tissue or surgical gloves.
- ▶ Check the product for loose or missing parts.
- ▶ Immediately put aside damaged or inoperative products and send them to Aesculap Technical Service, see Technical service.

3.10.2 Functional test

- ▶ Check that the product functions correctly.
- ▶ Check that all moving parts are working properly (e.g. hinges, locks/latches, sliding parts etc.).
- ▶ Check for compatibility with associated products.
- ▶ Immediately put aside inoperative products and send them to Aesculap Technical Service, see Technical service.

3.11 Packaging

- ▶ Appropriately protect products with fine working tips.
- ▶ Place the product in its holder or on a suitable tray. Ensure that sharp edges are covered.
- ▶ Package trays appropriately for the sterilization process (e.g. in Aesculap sterile containers).
- ▶ Ensure that the packaging provides sufficient protection against contamination of the product during storage.

3.12 Steam sterilization

Note

This product should be sterilised using Steam or Hot Air

- ▶ Check to ensure that the sterilizing agent will come into contact with all external and internal surfaces (e.g., by opening any valves and faucets).
- ▶ Validated sterilization process
 - Steam sterilization using fractional vacuum process
 - Steam sterilizer according to DIN EN 285 and validated according to DIN EN ISO 17665
 - Sterilization using fractionated vacuum process at 134 °C/holding time 5 min
- ▶ If several devices are sterilized at the same time in the same steam sterilizer: Ensure that the maximum permitted load according to the manufacturers' specifications is not exceeded.

3.13 Storage

- ▶ Store sterile products in germ-proof packaging, protected from dust, in a dry, dark, temperature-controlled area.

4. Technical service

⚠ CAUTION

Modifications carried out on medical technical equipment may result in loss of guarantee/warranty rights and forfeiture of applicable licenses.

- ▶ Do not modify the product.
- ▶ For service and repairs, please contact your national B. Braun/Aesculap agency.

Service addresses

Aesculap Technischer Service

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78532 Tuttlingen / Germany

Phone: +49 7461 95-1601

Fax: +49 7461 16-2887

E-Mail: ats@aesculap.de

Other service addresses can be obtained from the address indicated above.

5. Disposal

⚠ WARNING

Risk of infection due to contaminated products!

- ▶ Adhere to national regulations when disposing of or recycling the product, its components and its packaging.

⚠ WARNING

Risk of injury due to sharp-edged and/or pointed products!

- ▶ When disposing of or recycling the product, ensure that the packaging prevents injury by the product.

Note

The user institution is obliged to reprocess the product before its disposal, see Validated reprocessing procedure.

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