
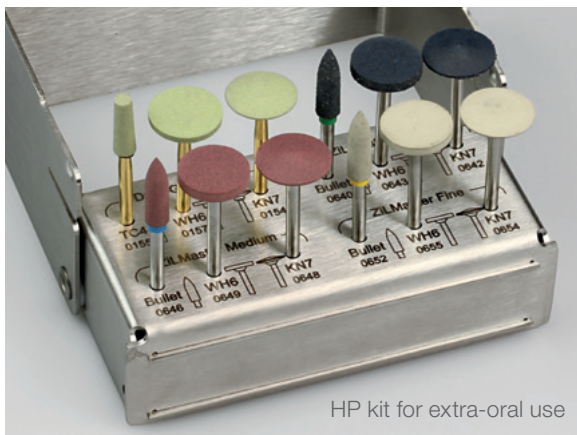


ZiLMaster Adjustment Kits

Time-efficient adjusting, finishing and polishing kits for zirconia and lithium disilicate

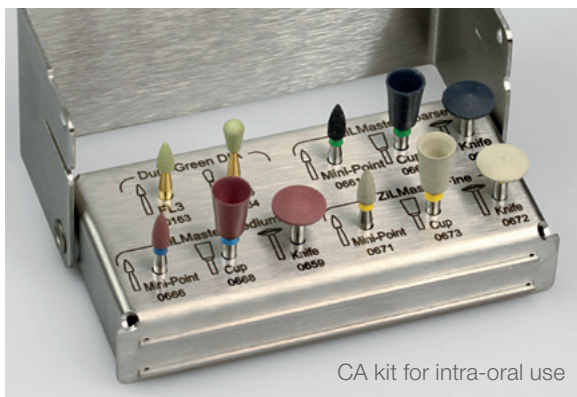


Smooth and homogeneous surfaces are essential to a harmonious integration of ceramic restorations into the natural dentition.

The instruments of the ZiLMaster Adjustment Kits have been specially designed for quick and safe contouring, finishing and polishing of monolithic CAD/CAM restorations.

Our kits include **Dura-Green Dia** diamond stones, developed specifically for monolithic materials, and **ZiLMaster** silicone polishers in three different grits, ideally tailored to their respective indications for use.

This well-matched step-by-step system is easy and time-efficient to use, from the correction of premature contacts to final super-polishing.

The instruments are available in two shank types, for contra-angle or straight handpieces, with a sterilisable bur block made of high-quality stainless steel.

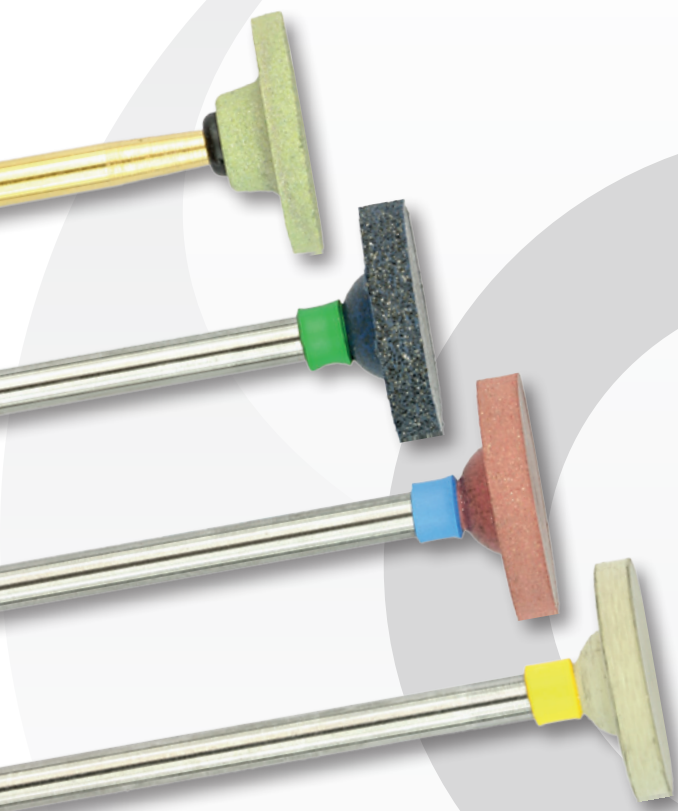

Dura-Green DIA



The diamond-impregnated microstructure of the Dura-Green DIA abrasives allows you to accurately and gently adjust contours or occlusal surfaces after correcting premature contacts, without any risk of chipping or micro-cracking.



Low contact pressure and smooth, low-vibration running properties make the instruments easy and comfortable to use.



The high edge stability ensures a long service life of the abrasives. The fine surface pattern created is an ideal basis for the subsequent polishing procedure.



ZiLMaster


The three different grits of these diamond impregnated silicone polishers ensure the optimum finishing and polishing of all-ceramic restorations.

Coarse (anthracite grey, green banded shank)



Suitable for abrasive prepolishing and finishing of adjusted areas or zones that still need minor material reduction

Medium (red brown, blue banded shank)

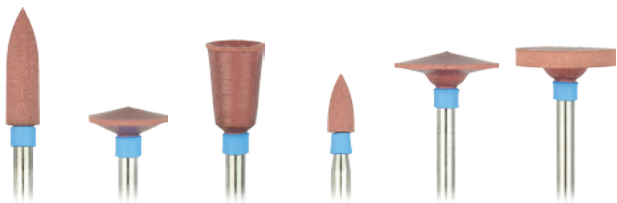


Recommended for contouring and smoothing anatomical details and preparing surfaces for super-polishing

Fine (light grey, yellow banded shank)



Designed for efficient super-polishing without the additional use of polishing pastes

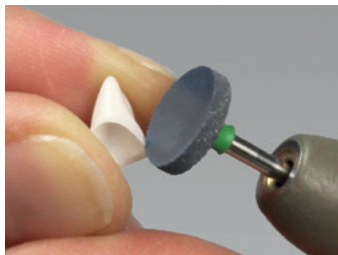

Minor corrections,
thinning margins and
adjusting contacts with
ZiLMaster Coarse



Smoothing surfaces,
marginal finishing and
polishing with ZiLMaster
Medium



Final super-polishing
with ZiLMaster Fine

highlights

Designed specifically for monolithic zirconia and lithium disilicate

Correction of premature contacts

Finishing adjusted contacts with low heat generation

Time-efficient high-gloss polishing

Surface improvement after insufficient glaze firing

Order information

Dura-Green DIA



| Shape | CY5 | FL3 | IC5 | KN7 | RD2 | RD3 | TC2 | TC4 | WH6 |
|-----------|------|------|------|------|------|------|------|------|------|
| HP 3 pcs. | | 0158 | 0162 | | | 0159 | 0161 | 0155 | |
| HP 2 pcs. | 0156 | | | 0154 | | | | | 0157 |
| CA 3 pcs. | | 0163 | | | | 0164 | 0165 | | |
| FG 3 pcs. | | 0168 | | | 0167 | | | | |

ZiLMaster



| | Shape | Bullet | Knife | Cup | Minipoint | KN7 | WH6 |
|--------|-----------|--------|-------|------|-----------|------|------|
| Coarse | HP 3 pcs. | 0640 | | | 0641 | 0642 | 0643 |
| | CA 3 pcs. | 0658 | 0662 | 0663 | 0661 | | |
| Medium | HP 3 pcs. | 0646 | | | 0647 | 0648 | 0649 |
| | CA 3 pcs. | 0665 | 0659 | 0668 | 0666 | | |
| Fine | HP 3 pcs. | 0652 | | | 0653 | 0654 | 0655 |
| | CA 3 pcs. | 0670 | 0672 | 0673 | 0671 | | |

Unmounted wheels are also available in all three grits.



ZiLMaster Adjustment Kit HP
shank
PN 0680



ZiLMaster Adjustment Kit CA
shank
PN 0681



Manufacturer
SHOFU INC.

EC REP SHOFU DENTAL GmbH

An der Pönt 70, 40885 Ratingen, Germany

11 Kamitakamatsu-cho,
Fukuine, Higashiyama-ku,
Kyoto 605-0983, Japan

www.shofu.de